INSTRUCTION MANUAL
Foot Massager | HMF-300

Mục lục

Safety cautions	01
Product structure	04
Product Features	05
Instruction before use	06
Operation guide	08
Removable cloth cover	09
Illustration	10
Maintenance and storage	11
Troubleshooting	12

Safety Instructions

Thank you for purchasing foot massager. Please read the Safety Precautions carefully before using the leg massager. The Safety Precautions are divided into 'warnings' and 'attention' in terms of the grade of danger. This book should be kept available for future reference after reading.

Important safety instructions

To ensure a safe and proper operation, please consult the instruction manual carefully before using this product.

- Use this product only as described in this manual.
 - Always unplug this product from the electrical outlet immediately after use and before cleaning.
 - Do not use this product if you have sensitive skin, neuropathy or poor circulation.
- Don't let children or people with mental disabilities use the heating function without proper supervision. Improper use may be dangerous.
- Do not operate this product under a blanket or when covered. Over-heating can occur and cause fire, electric shock or injury.
 - Do not cover this product in order to avoid electrical faults or damages due to excessively high motor temperature.
 - Stop the massage immediately if you feel any discomfort.
 - Never operate this product if it has been dropped, has a damaged cord or plug; experiences a malfunction, or if it has been dropped in water. Please contact an authorized service centre to do the necessary check and repair.
 - Do not use the components that are not recommended by Distributors.
 - Do not insert any objects into any openings in the product.
 - If you want to cut off the power supply, please make sure you have turned off all the operating controls, as well as remove the power plug from the socket.
 - Please do not use this product as a medical device. If you need a professional diagnosis or treatment, please consult a professional medical person.

Safety Instructions

Operating Environment Precautions

Do not use this massager in humid or extremely dusty environments to avoid malfunction.

- Do not use this massager at a temperature above 40 degree centigrade.
- Do not use this massager outdoors.
- Avoid placing or operating this massager in direct sunlight or near a heating appliance, stove

User precautions

Make sure that this product is far away from any hot surfaces.

- Please consult your doctor before use for those who have the following diseases:
 - Those who have heart disease, malignant tumor, back pain, spinal scoliosis, hypertension, osteoporosis, etc.
 - Those who are pregnant or having a menstrual period.
 - Those who have a physical discomfort, fever, a swelling or inflammation, bruises, allergies, etc.
 - Those who are warned by the doctor not to do any massage to avoid injury or physical discomforts such as blood clots (thrombosis), severe aneurysm, acute venous tumor, all kinds of skin inflammation, infectious skin disease (including subcutaneous tissue inflammation), etc.
 - Warnings: Do not use the massager together with the following medical instruments.
 - 1) Those who are using a pacemaker or other implantable medical equipment.
 - 2) Those who are using an artificial heart.
 - 3) Those who are wearing an electrocardiograph for similar medical equipment.Otherwise it will cause malfunction for those medical equipment, putting the user's life in danger.

Safety Instructions


Precautions for power supply

- Do not turn on the switch when your hands are wet.
 - Do not immerse in or spray water or liquid on the control power switch.
- e) Cord/Plug Precautions
- Do not connect or remove the plug with wet hands. Otherwise you will suffer an electric shock or a product malfunction.
 - Do not pull the wires when you pull out the power plug. Please hold the plug itself.
 - As the default operating time is 15 minutes, all the functions will stop working after 15 minutes. Please unplug this product from the electrical outlet immediately after use.
 - Keep the cord away from hot surfaces.
 - Do not place the power plug under the product or any heavy items.
 - Do not wind the power plug around the product lest it should be damaged and cause an electric shock or fire.
 - Do not operate this product if it has a damaged cord or plug.
 - Do not use this product if the power socket is loose.
 - Do not carry this product by its supply cord or use the cord as a handle.
 - When the power plug is damaged, please call the service centre or an authorized technical person for the necessary check and maintenance.

Attention

- Do not fall asleep while massager is on.
- Do not use this product while you are drunk or physically uncomfortable.
- Don't stand on the product or put a heavy weight on it.
- Do not use this product for more than 30 minutes.
- Do not use this product if the cloth cover is worn out.
- Do not use this product within an hour before or after a meal.
- When you are using the rollover function, when the flip angle has reached its maximum, please do not continue to flip hard, or it may damage the product.
- When you are using the product at a low room temperature, do not suddenly use the heat function. You can adjust the room temperature to normal room temperature, before starting the heat function.
- When moving this product from a low temperature storage place to a warmer setting, it may cause malfunction due to water vapor gathered in the components. Set the product in a normal room temperature for one hour before using it.
- Do not use this product continuously if you find any components damaged or the product in an abnormal condition.
- When current interruption occurs during your massage, try to avoid further damage due to sudden power restoration.
- Stop using the product if there is thunder and lightning.

Product structure

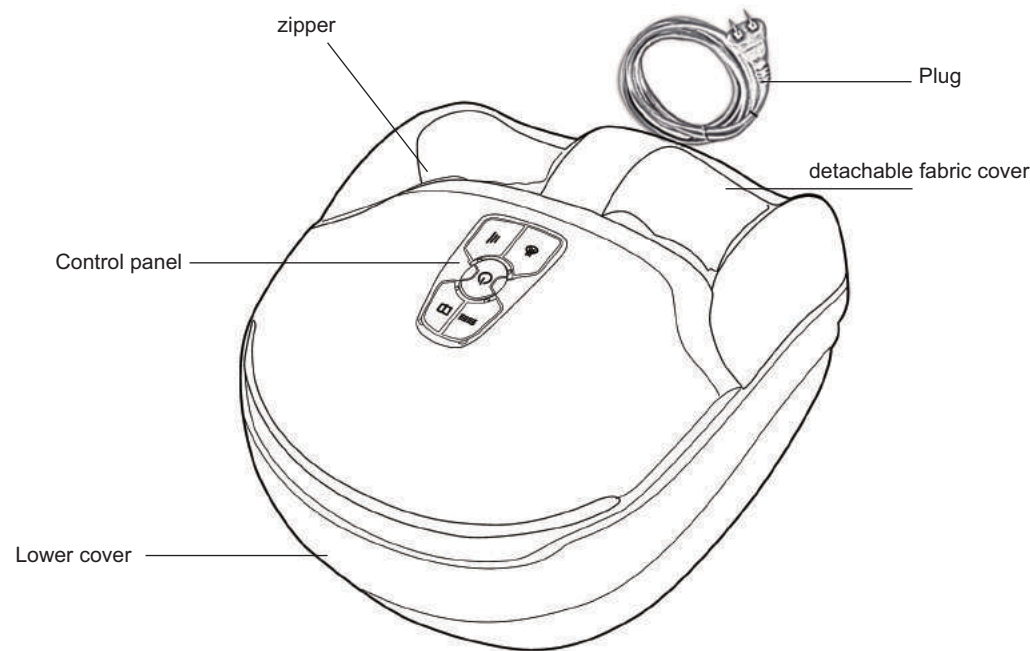Product Features

New innovative design

new design, car shape, stylish and portable .

Full coverage of massage

Combination variety of functions of shiatsu,kneading,air compression,heating, vibration to provide deep kneading shiatsu on nodes and stimulate all points on foot reflexology zone.

New massage modes

Vibration on feet and legs to alleviate muscle pain.

Detachable cloth cover

Zip-out cloth cover on feet for easy cleaning and changeable.

Safety

In the process of massage, in case power is off suddenly, the airbags will be deflated automatically without holding feet tightly.


Product Features

Relieve foot, leg pressure, improve muscle flexibility, beautify the leg line, promote blood circulation, relieve tension pain after exercise, stimulate important points, improve overall health, relieve the nervous system, reduce fatigue and improve sleep.

Instruction before use

1.How to place the machine

please don't expose the machine to the direct sunlight put the carpet under the machine to lower working noise.

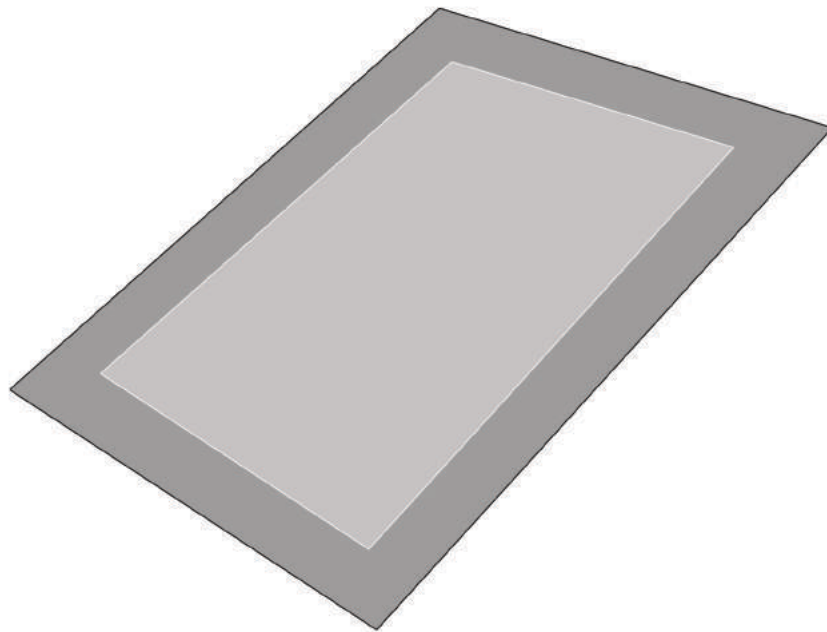How to use it

make sure the surface of the wire is okay, Do not insert or remove the plug with a wet hand

This product uses only the DC110V-240V, 50/60Hz power supply. Make sure your indoor voltage fits this specification parameter.

- Connect the DC adaptor to the socket.
- Turn on the power switch at the back of product to "I", when the red light is on, the yellow light on control panel blinks, the product is in standby mode.
- Press on/off button, start the foot massage.
- Chose different program to proceed massage.


Product Operation Guide

Product Operation Guide:

Auto program modes: Place your foots in the product, before pressing the switch key on screen. When the blue light is on is mode 1, Red light is on is mode 2, Green lights on is for mode 3. program.

Vibration: Press on/off, It will control vibration function. Lights on, vibration function on, Lights off, vibration function is off. Under stand by mode, press vibration for 3 seconds, it can solely control vibration. 3 modes of vibration adjustment, Blue lights on for mode 1. Red lights on for mode 2. Green lights on for mode 3. mode 3.

Heat : Press the button, Red lights on for heating function. Lights off for close heating function.

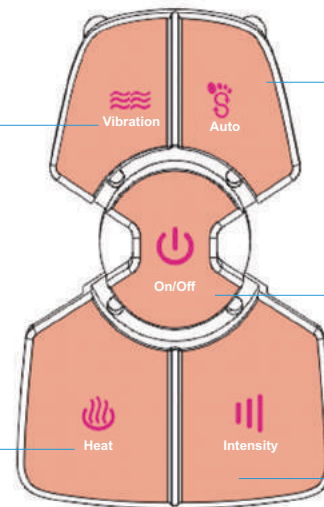On/off : Start the foot massage, turn off the foot massage. When you start the massage, the lights will flash. Press it, green lights on, the massage function is turned on, and start vibration, automatic program, heat function is on. Press it again for close massage function.


Intensity: 3 levels of air-compression intensity adjustment, green lights on for low level. Blue lights on for media level. Red lights on for strong level. Light off for close function.

Note: 1, when you turn the vibration function keys to adjust the intensity of the three levels of vibration.
2, under automatic massage function, intensity key is for adjust the intensity of air-compression.

Removable cloth cover


Take off the cloth cover

- 1.unzipper it
- 2.take off the cloth cover.

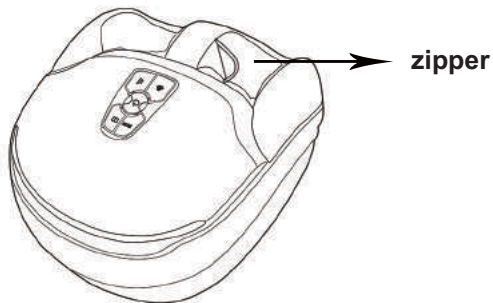
Load in cloth cover

- 1.Attached the sticker
- 2.Put on the zipper.

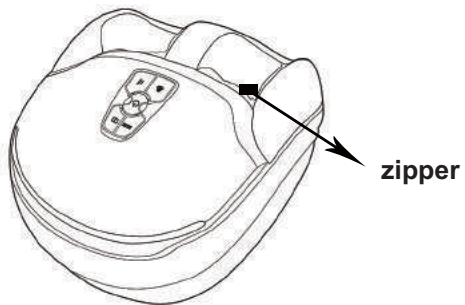Illustration

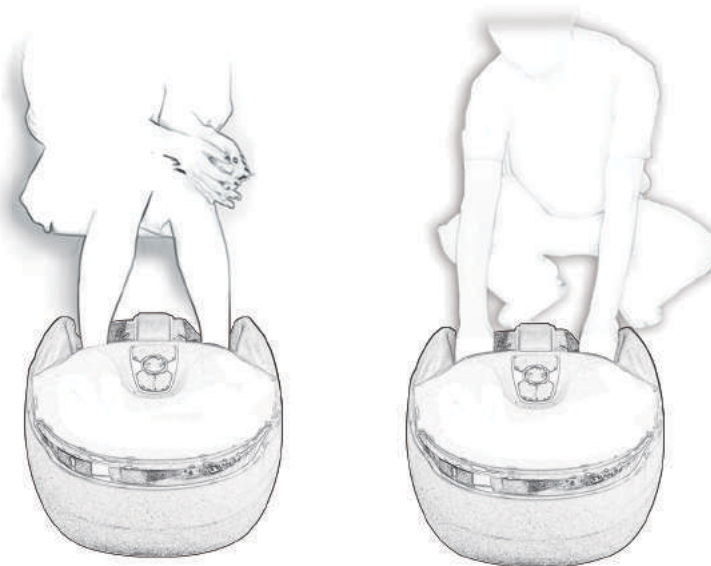Do not put your hands inside the cloth cover.

Maintenance and storage

Product Maintenance

Storage

- Please store it in the packing carton, and put in a safe, dry, well-ventilated, and shaded place.
- Do not twist the power plug.
- Do not hang the power wire and pull the product.
- Avoid sharp items from destroying the surface of the product.
- Please don't store the machine under direct sunlight or in places of high temperature.

Clean

- Before cleaning, make sure that the power cord is unplugged from the electrical outlet and that the machine is completely cool.
- Use a damp cloth to clean the product. Ensure the product is completely dry before reuse.
- The detachable fabric cloth can be washed

Attention

Make sure not to use benzene or thinner and other corrosive detergents to clean the product.

Maintenance

Keep the product clean and store in a dry environment.
When the product is not in use for a long time, please put a dustproof cover on the

Troubleshooting

If you meet the following problems, please referring to the following solutions. If the product is still malfunctioning or doesn't work, please send the product to our service centre or an authorized technician for the necessary checks and repair.

Conditions	Possible reason	Solutions
It doesn't work	1. Plug isn't connected correctly to the socket. 2. The power switch at the back is off. 3. Switch key on control panel isn't on. 4. The fuse blows out.	1. Push the plug firmly into the socket. 2. Switch on the power switch at the back of the product. 3. Switch on the on/off key on the panel. 4. Insert a tooth pick in the small hole of "fuse" under the switch power, before taking out the fuse.
Stop working suddenly	1. The power switch is off. 2. Massage time is over 15 minutes. 3. If it is working for a long time, the motor will be protected against overheating. 4. If you put heavy pressure on the rollers and the rollers do not work well.	1. Check that the power plug is connected well. 2. Restart the power supply switch. 3. After the product is cooled down, please restart the ON/OFF power. 4. Switch off the power supply, and turn it on again.

Product Specifications

Product name	FOOT MASSAGER
Model number of the product	HMF-300
Voltage	AC 220V,50/60Hz
Power consumption	65W
Timer	15minutes

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