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
## GHẾ MASSAGE
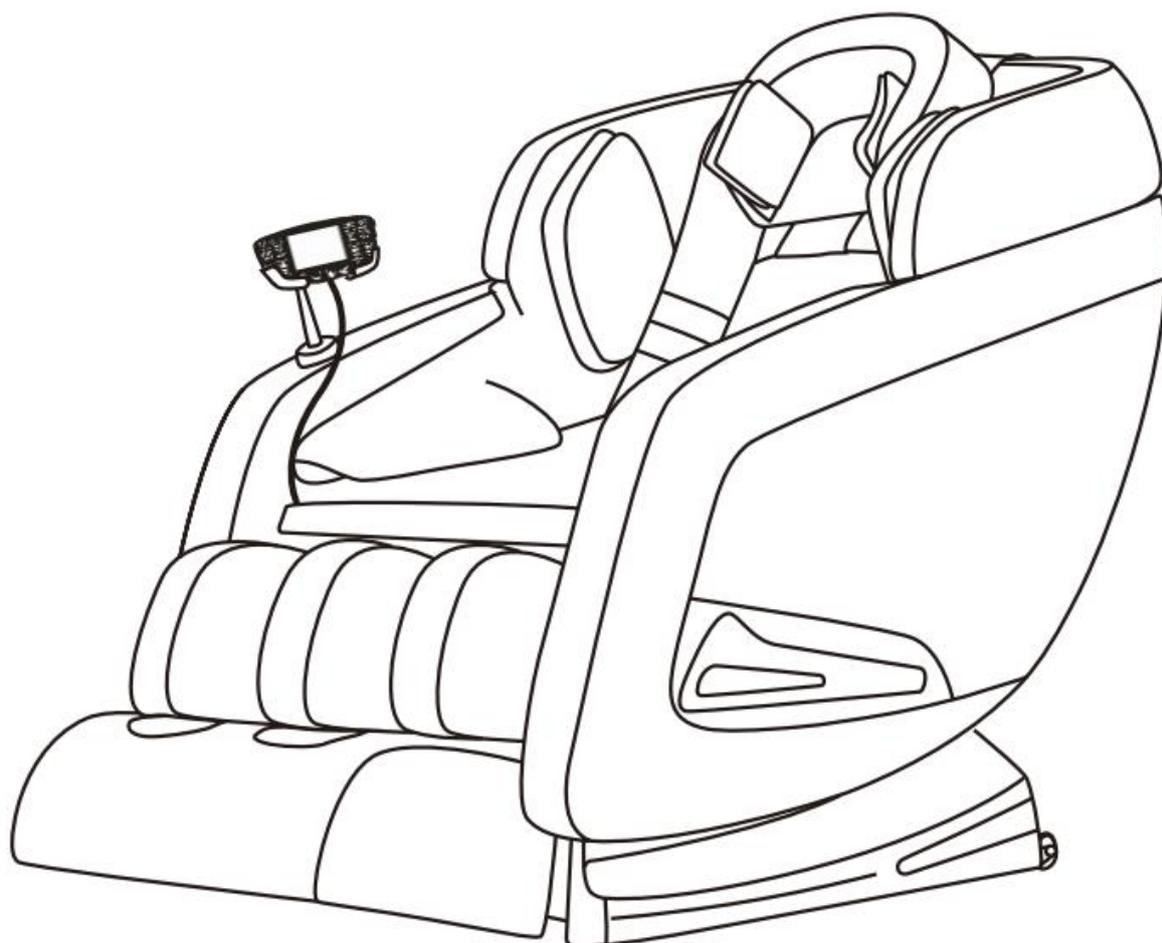
# INSTRUCTION MANUAL

## MASSAGE CHAIR


**Model: MK-331**

**Read this instruction manual carefully before using!**

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## GENERAL DESCRIPTION
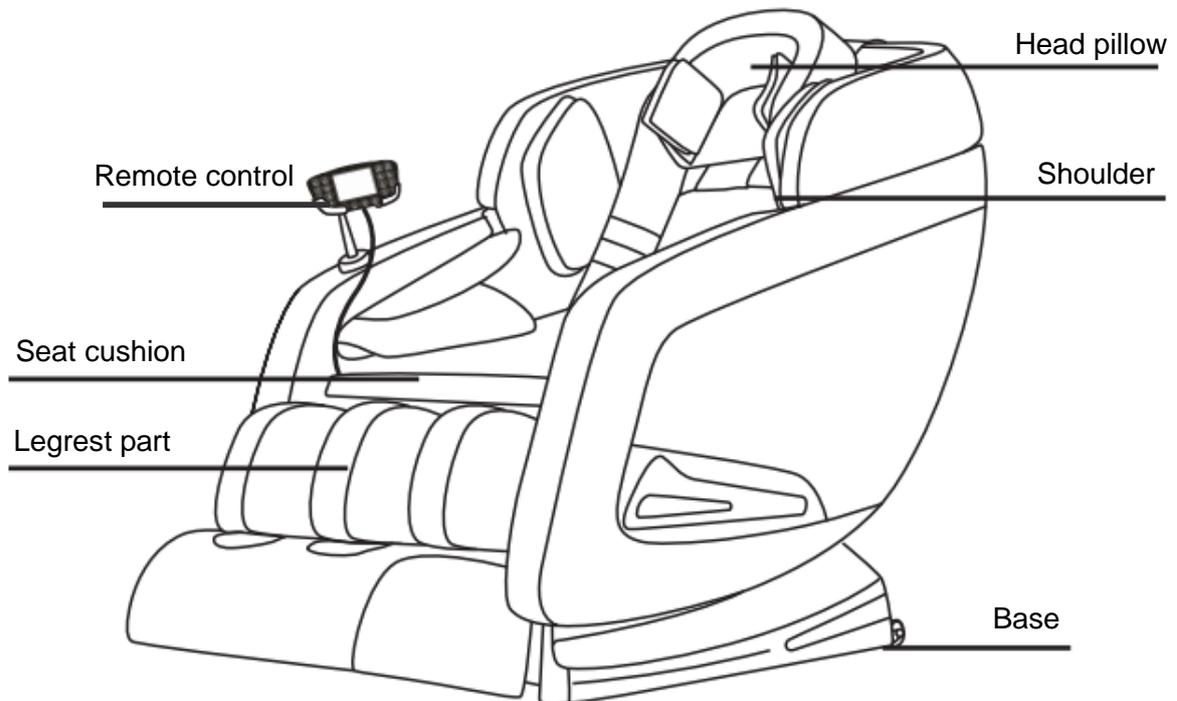
## IMPORTANT SAFETY INSTRUCTIONS

Please read all instructions before use. The massage chair must use a socket that matches the plug. Incorrect grounding will result in electric shock. Please do not change the plug and socket without permission. In the use of electrical products, in order to reduce the risk of electric shock, burns, or personnel injuries, basic safety precautions shall be observed, including the following:

- Be sure to use the socket that matches the plug. Make sure the plug is completely plugged into the socket to reduce the risk of short circuits or fire.
- Do not store or use the product outdoors or in direct sunlight or under high humidity.
- If there is a sudden power failure during the operation, please switch the switch to the opposition immediately and unplug the plug so as not to damage the massage chair.
- Do not operate massage chairs if the wires are wet or damaged.
- Please do not use wet hands to insert and remove the plug so as to avoid accidents.
- When a child, a handicapped person or an invalid is near a massage chair or using a massage chair, someone must be there.
- This product can only be used for purposes described in this manual.
- Do not enter the sleeping state during the massage.
- Always keep the massage chair clean and avoid inserting other objects into the gap of the massage chair.
- Please check the seat back, lint and so on before each use. If abnormal, please do not use and contact the customer service as soon as possible.
- To avoid discomfort or injury caused by the use of this massage chair, please consult a

doctor before using this product.

- People who cannot use the massage chair: A person suffering from malignancy, suffer from heart disease; pregnant woman, or period person; people whos being treated; people with a back disorder or abnormal or distorted cervical spine.
- If you feel pain during use, discontinue use immediately. If you are injured during use, seek medical attention immediately.
- This product is only for family comfort massage enjoyment, not professional massage medical instruments, therefore, can not replace professional treatment.
- It is recommended to use a massage chair for 30 mins, massage the same position of not more than 10 minutes, otherwise excessive massage will lead to counterproductive effect, and long time use can cause the massage chair overheating and shortening product life.
- Do not use this product in gas or oxygen environment.
- Children under 18 years of age are not eligible for this product

**Warning:** This product is intended for those who are physically, mentally, mentally handicapped, or inexperienced, unless they have received the guardianship or instruction of the person in charge of the product.

Children should be supervised to make sure they do not play with this massage chair.

If the cord is damaged, it must be replaced by the manufacturer or its maintenance department or similar professional personnel to avoid danger.


Equipment surface heating, heat insensitive people must pay attention to the use.

## PRODUCT SPECIFICATIONS

Product name	Massage chair	Product size	122 x 74 x 110 cm
Model	MK-331	Net weight	60 Kgs
Rated voltage	AC220V/ 50Hz	Gross weight	70 Kgs
Rated power	90 W	Packing size	119 x 77.5 x 109 cm


## PRODUCT INSTALLATION

- Open the packing box and take out the parts as below:

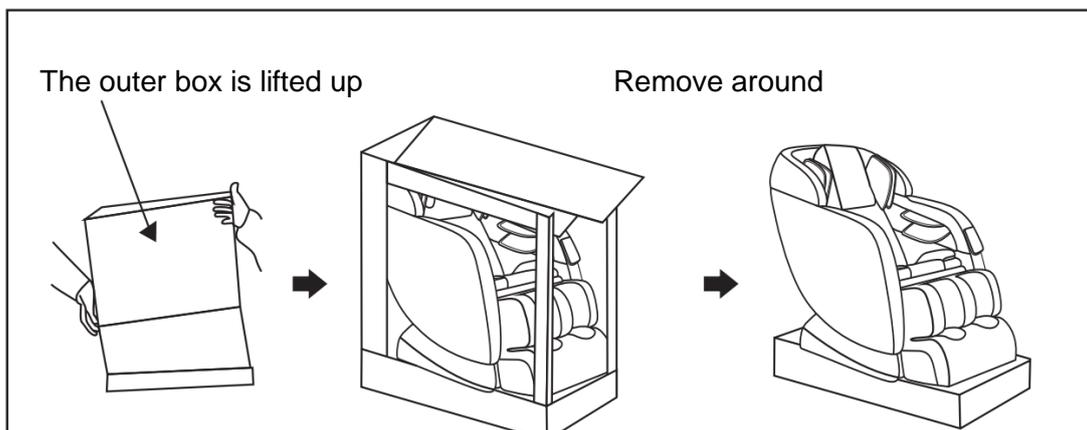The massage chair is fully packed and free from technical installation.

- Place the massage chair:

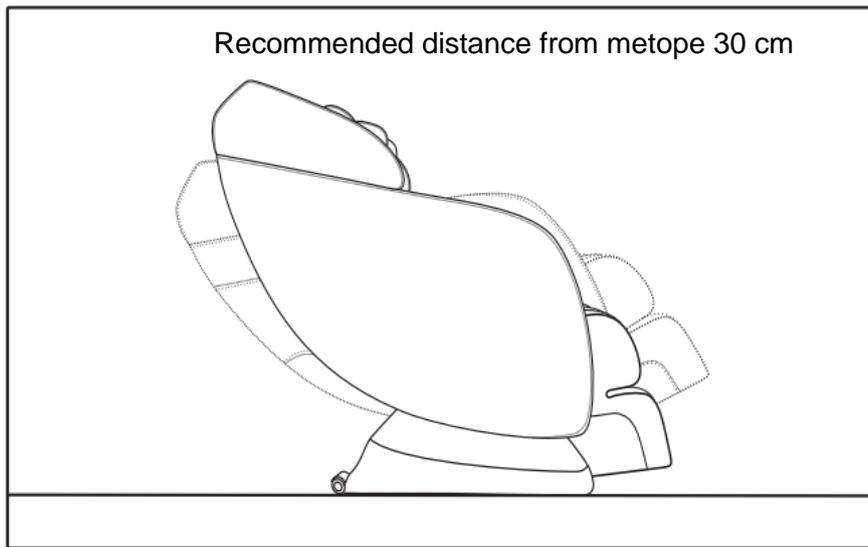
Lift the whole massage chair out of the bottom box and place in the right place.

When placing the massage chair, be sure to leave 30 cm space on the back of the chair. Leave 55 cm space in front of the chair so that the legs and backrest of the chair can be freely dumped.

\*Do not put the power cord in the electric control box of the massage chair.

\*In order to prevent ground damage and reduce noise, please put the cloth under the massage chair and so on.

- Moving the massage chair:

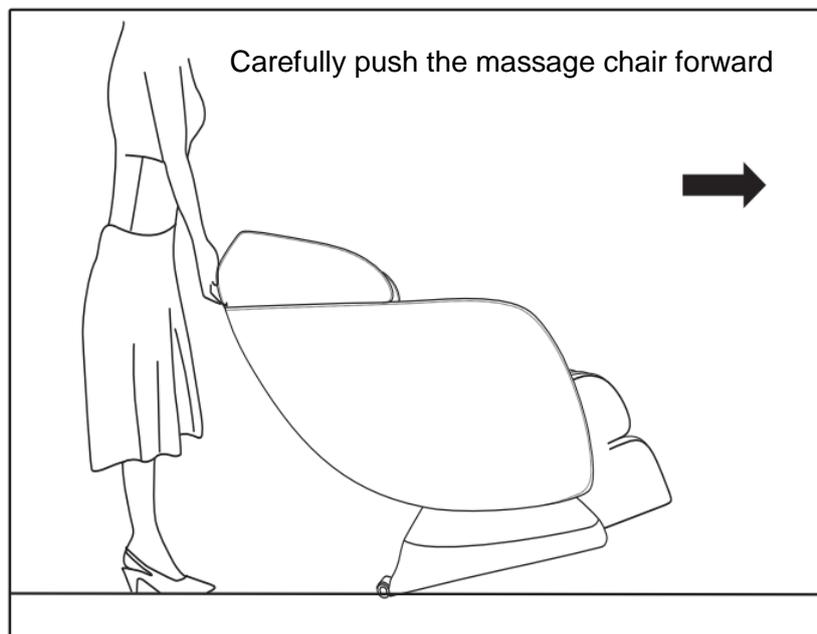As shown in the first with both hands on the chair backrest, dumped back 45 degrees, be careful pushing backward or forwards casters massage chair, until the right down position.


When someone is lying on the massage chair, do not move the massage chair. It may damage the floor when moving. It is recommend to lower the mattress on this product.

In order to avoid falling, please firmly hold the massage chair when you move it.

# OPERATING INSTRUCTION

## Control panel

### Control panel function zone: Controller front

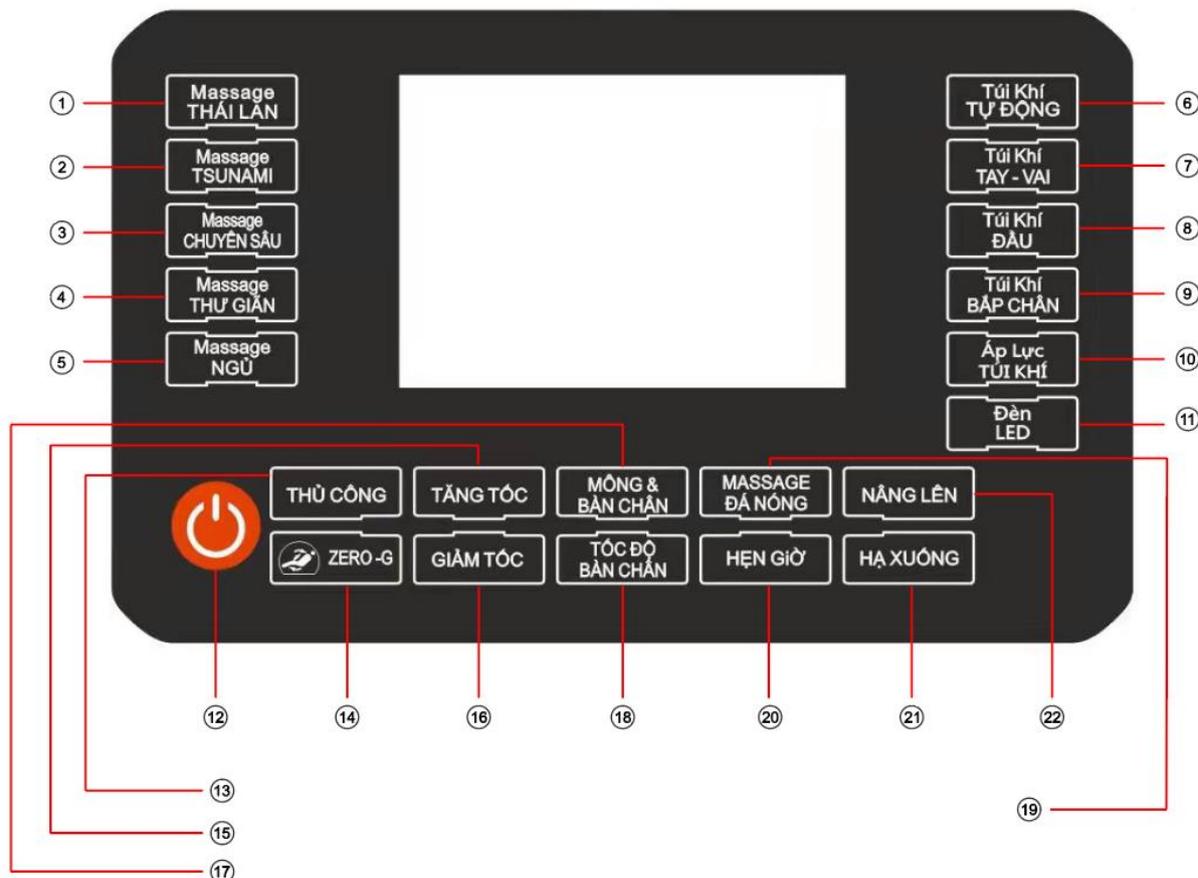- |                        |                        |                   |
|------------------------|------------------------|-------------------|
| 1. Thai massage        | 9. Leg airbag          | 17. Buttocks feet |
| 2. Tsunami massage     | 10. Pressure intensity | 18. Feet speed    |
| 3. Strong              | 11. Led light          | 19. Heating       |
| 4. Relax               | 12. Power              | 20. Timing        |
| 5. Sleep               | 13. Manual mode        | 21. Drop          |
| 6. Pressure auto       | 14. Zero gravity mode  | 22. Rise          |
| 7. Shoulder-Arm airbag | 15. Speed +            |                   |
| 8. Seat head airbag    | 16. Speed -            |                   |

### Control panel button function description

#### Power/On/Off

Press to turn the massage chair on or off.

Press the power button and select either automatic massage or manual massage option. The base time is set to 15 minutes.

**Automatic mode program selection: Thai massage (Kneading), Tsunami, Strong, Relax, and Sleep**

**Pressure auto:** Press to enter to air pressure auto mode.

**Shoulder Arm Airbag:** Press to turn the shoulder and arm airbag on or off.

**Seal Head Airbag:** Press to turn the seal cushion head airbag on or off.

**Leg Airbag:** Press to turn the leg airbag on or off.

**Pressure intensity:** Air strong regulate.

**Zero gravity mode:** Press to move the chair into the zero gravity position (3 steps).

**Foot speed:** Press to turn the foot speed on.

**Speed+/-:** Press this button when you want to increase or decrease the massage speed.

**Timing:** Use this button to change the massage time anywhere between 5 and 30 minutes. The default time is 30 mins.

**Led light:** Press this button to turn on the light on side.

**Heating:** Press to turn the heating mode on or off.

**Rise:** Long press, calf rise, backrest down.


**Drop:** Long press, lower leg, backrest rises.

## TROUBLESHOOTING PROBLEMS

In case of an abnormal situation, please stop using and unplug the power supply to prevent the accident. Please contact the service center if the following process is still invalid.

Possible anomalies	Processing method
The sound produced during operation.	The sound of the structured design does not affect the machine
Plug in the power supply, but the massage chair is inoperable.	Check whether the power plug is plugged in or whether it is energized. Check the power on the massage chair to make sure the key is in the "open" position.
The use of the process, the controller can operate, but the massage chair can not be operated.	Please turn off the power supply and stop using for 30 to 50 minutes, so that the inner heat of the machine will be cooled.
In use, the controller is inoperable and the massage chair stop operating.	Auto timer function, after starting, if the set time has arrived, please press the button again.

## CLEANING AND MAINTENANCE

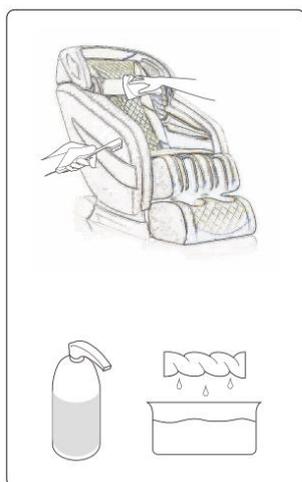- \* With a dry cloth to wipe the sweat.
- \* Use low concentration of neutral detergent when cleaning stains.
- \* Soak the soft cloth, wring it thoroughly, and gently wipe the massage chair.


\* The power plug must be removed when cleaning the massage chair.

\* In addition, do not touch the power plug with wet hands.

\* Please don't clean the massage chair with water.

\* Do not use a blower to dry the sweat on the massage chair.



\* Do not use gasoline, thinners, alkaline washing, bleach,... to clean.
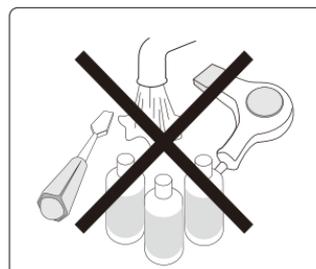
# INSTRUCTION OF SCANNING QR CODE

## Method 1: Scanning QR code by Zalo (Vietnam app)

- Open Zalo.
- Click ( + ) and select Scan QR code

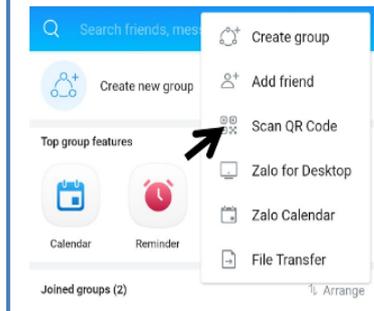
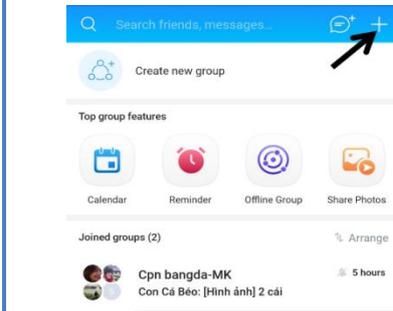
- Move the camera to QR code
- Click "Yes" to redirect to the warranty registration page

## Method 2: Scan QR code on iPhone with iOS11.

- Activate scanning QR code by camera: Go to Settings > Camera > Scan QR Code).

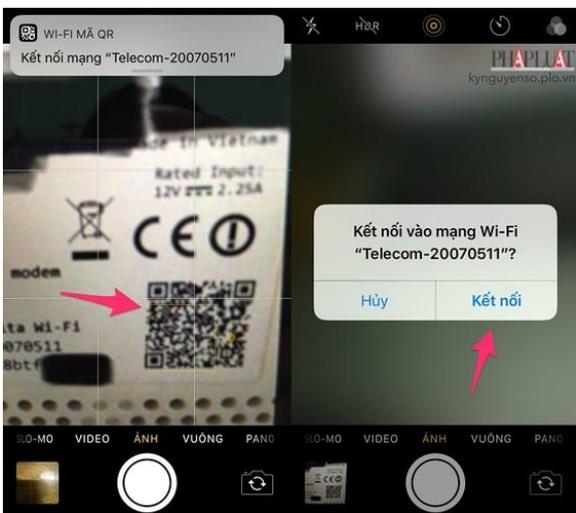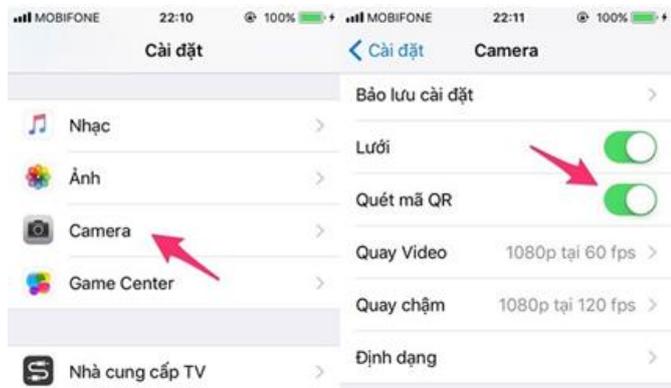- After completed install, open camera and move the camera to QR code, then you will receive a notification, please click "Connect" and follow the instructions.


### Method 3: Scan QR code on Android with Google Assistant

Instructions for using Google Assistant on Android.

Hold the Home button to open the Google Assistant application and follow the illustration

After clicking on the icon above the application to open the camera, move it to QR code and click on the link in the image

Click the "Start" and follow the instructions
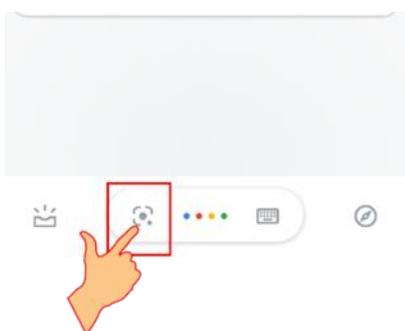
### CUSTOMER SERVICE CENTER

South: 12/9 Tay Quarter, Highway 13, Vinh Phu, Thuan An, Binh Duong Province.

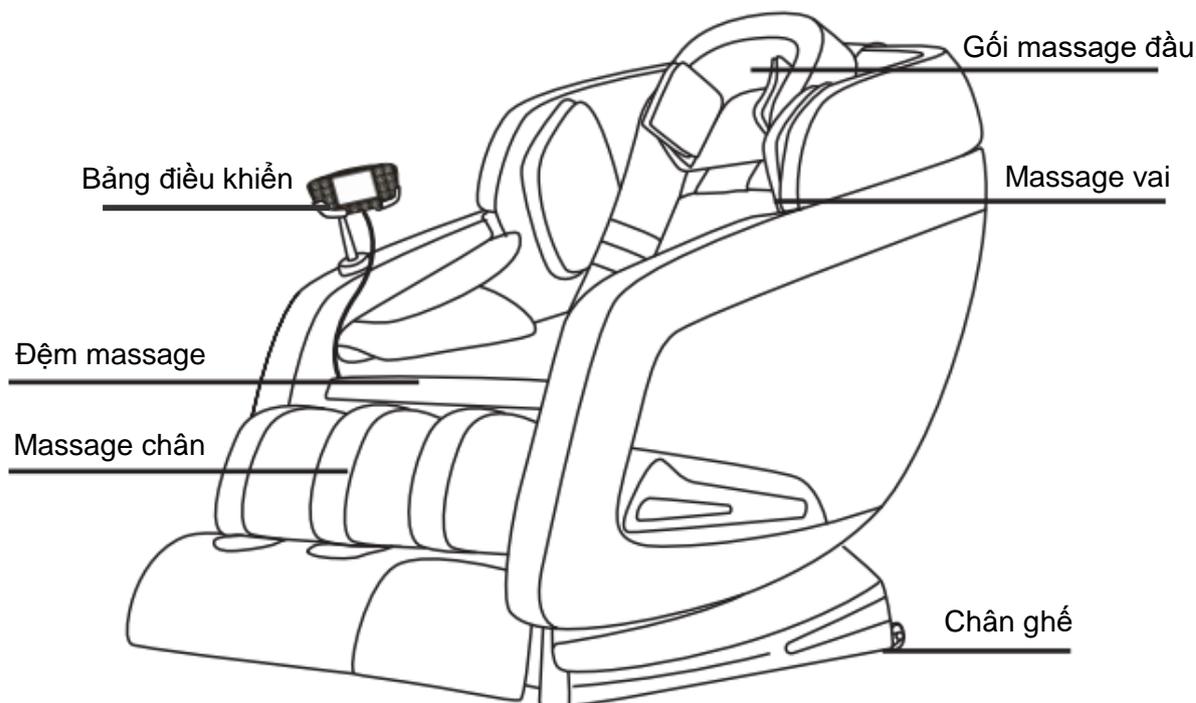
North: Xam Thi Village, Hong Van Commune, Thuong Tin District, Hanoi

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