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
SÚNG MASSAGE

MODEL: MK377

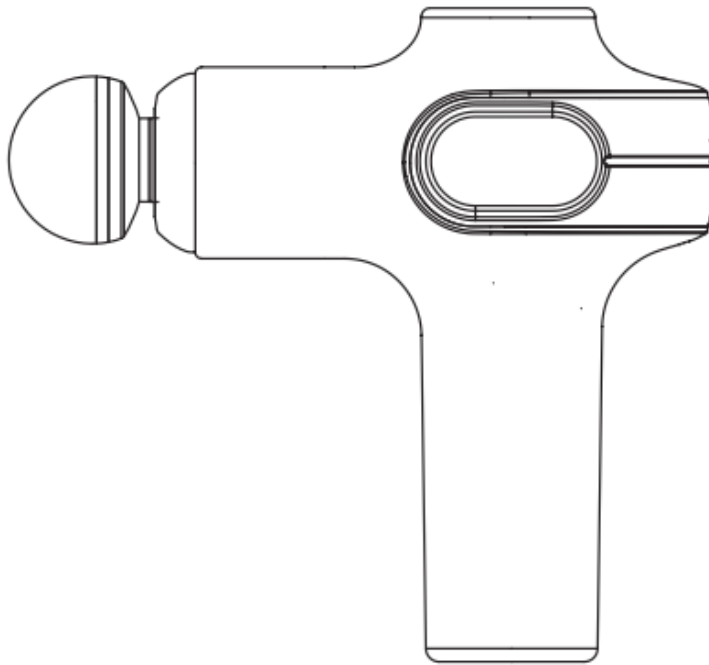INSTRUCTION MANUAL

MASSAGE GUN

MODEL: MK377

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Read this instruction manual carefully before using

*****Technical specifications :**

Product Name	Massage gun
Model	MK377
Voltage	DC7.4V
Rated power	18W
Dimensions	190MM*165MM*60MM

*****Using Instructions**

Please insert the massage head into the front hole of the machine before use.

- ① Power on: Press and hold the power button, the indicator light will light up, and the power on is complete.
- ② Gear shift: After turning on the machine, briefly press the power button to increase the gear and gradually increase the machine's force from 1st gear to higher gears. Press to return to first gear after reaching the highest gear.
- ③ Shutdown: Press and hold the power button, the indicator light will turn off, and the machine will shut down.

*****Battery Instructions**

The battery needs to be charged for about 3 hours to fully charge when fully discharged, and the single endurance of the machine depends on the pressure level during use. When charging the battery, please use the specially equipped charging cable of the machine and match it with a 5V adapter that meets national standards and does not exceed 15W to charge the device.

When the machine is charging, the light flashes, and the light remains on when fully charged. When the battery is low, the light flashes five times and then goes out.

*****Attentions**

- The deep penetration of the device can relieve pain, bursitis, muscle spasm, neuralgia and other skeletal muscle diseases caused by joint sprain in a short time. Almost all soft tissue diseases are related to trauma, muscle overuse or abnormal posture. Most of these diseases involve facial scar tissue, which can be treated with the device. When using the equipment, excessive pressure and long-term physical therapy should be avoided.
- Sensitive areas such as head, face and cervical spine should also be avoided. If there is inflammation, swelling or aggravation of pain, the use of the equipment should be stopped immediately. The current research is about the effect of vibration training on specific medical diseases. This is likely to reduce the types of contraindications mentioned above. Practical experience shows that in many cases, the

combination of vibration training with physical therapy is also a physical therapy suggestion. This must be done on the advice of, or in the company of, a doctor, specialist or physiotherapist.

- In order to reduce the risk of failure, electric shock, injury, fire, death and equipment damage, please read all instructions and warning tips in this manual before using this product.
- This product is not used for medical purposes, only provides personal massage function
- Do not allow children or people with no autonomy to operate and use this product, please place this product out of children's reach;
- It is forbidden to use it on the surface of non-dry, clean skin, it can be used through clothes, press and move gently on the skin, about 10 seconds for each part;
- It is forbidden to use the product on the top of the head or any part of the body that is hard or has bones; if there is no pain or discomfort, use the product only on the soft tissues of the body;
- It is forbidden to use it when there is bruise, pain or discomfort in the massage area;
- Do not stand or sit directly on the product, it may cause a fall accident; Excessive use will cause the product to overheat and shorten the product life. When using repeatedly, let the product cool down for a period of time before using.
- Do not use it in a humid environment such as a bathroom, and do not put the product in water or other liquids to avoid short circuit or electric shock;
- Do not place the product near the stove or direct sunlight;
- Do not use during lightning strikes;
- Do not insert or drop objects into any openings of the product;
- Do not cover the massager when using it, so as to avoid malfunction due to the increase of motor temperature.
- Do not allow fingers, hair or other parts of the body to be close to the shaft and back of the vibrating head, as they may be caught.
- The battery must be removed from the appliance before it is discarded;
- When taking out the battery, the appliance must be powered off;

*****Maintenance**

1. Wipe the surface of the equipment with a slightly damp towel, and then dry it with a soft cloth.
2. When storing for a long time or traveling, please put the product in the original packing box of the equipment and keep it dry.

*****Warning**

- Please use our specially equipped USB cable with 5V power adapter to charge the product.
- Do not use a charger (adapter) with a damaged casing or an incorrect voltage to charge the unit.
- Do not touch the charging cable and power adapter with wet hands, or pull the charging cable to remove the charging cable or charger (adapter).

- Do not touch the product or charger (adapter) with wet hands to avoid product short-circuit, malfunction or electric shock. Please stop using the adapter when it is exposed to rain, liquid or severely damp.
- The following persons: implantable cardiac pacemaker, patients with heart disease, strain back disease, pregnant women or shortly after delivery, fracture, osteoporosis, joint disease, patients with malignant tumor, skin trauma or eczema, etc. are prohibited to use. If you are not fit, please consult your doctor before using
- Do not expose the battery to high temperature or around the heating equipment, such as sunshine, heat, etc. Overheated batteries of heaters, microwave ovens, ovens or water heaters may cause explosion.
- Do not disassemble or refit batteries, insert foreign matters, or immerse them in water or other liquids, so as to avoid leakage, overheat, fire or explosion of batteries.
- Before abandoning appliances, batteries must be removed from appliances.
- When removing batteries, appliances must be powered off.
- Batteries should be disposed of safely

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