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
MODEL: MK390

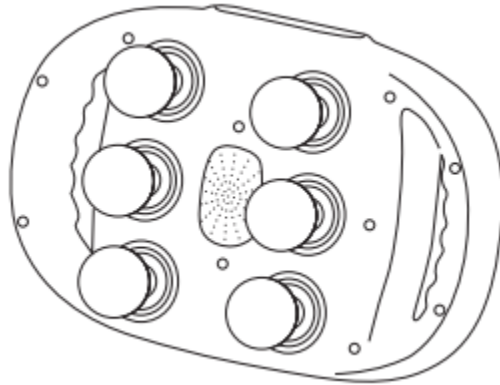USER MANUAL

SIX HEADED MASSAGER

MODEL: MK390

Read this instruction manual carefully before using

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***Technical specifications :

Product Name	Six headed massager
Model	MK390
Voltage/ Rated power	DC 5V/ 15W
Dimensions	27 x 15.5 x 11.5 cm

***Functional indication

① Type-C charging port


② Power on and mode switch keys

1. Long press power on and short press mode switching.
2. In boot mode, display the digital display for 10 minutes.
3. Short press to switch the heating light on or off, displaying P1 or P2
4. Continue short pressing: Mode P3 P4 P5 P6 P7 P8 P9

③ Digital display (white light) with two digits

④ Intensity increase button, L1-L9, after 2 seconds, the machine will return to the countdown interface

⑤ Intensity decrease button, L9-L1, after 2 seconds, the machine will return to the countdown interface

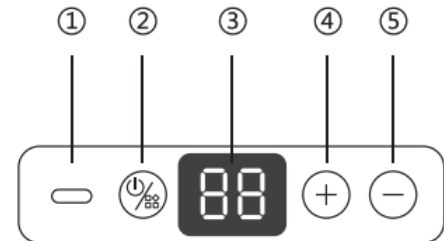When charging, the white number flashes, and when fully charged, the white number stays on for a long time

***Maintenance

1. When not in use for a long time, maintain a certain amount of electricity, turn off the power switch, and properly store the product.

It can be placed back in the original packaging box in a safe, dry, ventilated, and low-temperature place.

2. When not in use for a long time, it is recommended to charge the product every six months to prevent damage to the battery due to low power.

3. For daily maintenance, please use neutral cleaning agents and soft cloths to wipe. Do not use engine oil, organic solvents, or other chemical agents for cleaning, and do not bleach in water.

4. Before using or cleaning the product, please turn off the power switch.

5. During transportation, this product should be securely secured to prevent severe vibration, collision, rolling, and falling.

***** Notice**

1. It is prohibited to use a pacemaker, facial cleanser, or central venous catheter inserted through peripheral veins.
2. All kinds of acute and chronic infectious diseases, acute high fever, acute osteomyelitis, tuberculous arthritis, acute abdomen or infectious skin diseases, all kinds of skin ulcers, sores, burns in the acute stage, severe diabetes patients or recent fractures or wounds are prohibited.
3. It is prohibited for children under 6 years old to use it. Adolescents over 6 years old but under 18 years old should use it under adult supervision. During pregnancy and breastfeeding, women should avoid using the product on their body or under the guidance of a therapist.
4. Hemorrhagic diseases such as thrombocytopenia are prohibited from use. It is prohibited to replace or supplement medical treatment with this device. If you have pain, are receiving medical treatment, or use other devices, please consult your physician.
5. Before use, please confirm whether the voltage of the product is consistent with the ambient voltage, and whether the power plug is suitable for the socket
6. It is prohibited to immerse the product in water. If it accidentally enters the water, please send it for inspection in a timely manner.
7. It is prohibited to bring the product close to a fire source or put it into a fire source.
8. It is prohibited to use non original chargers, power cords, or plugs.
9. It is prohibited to intentionally drop this product.
10. It is prohibited to disassemble screws or attempt to disassemble this product.
11. If the product is damaged or has abnormal sounds, please do not use this product.
12. Do not use this product when your hands are wet or near water sources. For example: kitchen, bathroom.
13. When using, do not press or block the muzzle of the gun. Please do not apply heavy pressure to the fascia gun.
14. When the power cord or plug of the product is damaged, incomplete, or overheated, please do not disassemble or repair it without authorization. Please promptly contact the dealer or designated after-sales service for repair.
15. The appliance has a heating surface. Attention must be paid when using it for people who are not sensitive to heat.

***** Frequently asked questions**

Unable to start or unable to start due to low temperature

Reason 1: low temperature

(Prolonged exposure to environments below 0 degrees Celsius can lead to decreased battery activity; failure to start properly, resulting in the entire machine appearing to be able to start normally but unable

to operate normally)

Solution: Let it stand still for 2 hours in an environment above 0 degrees Celsius, or use a hair dryer to quickly increase the body temperature.

Reason 2: Insufficient battery power

(When the battery is in a low battery state, if it is not charged for a long time, the battery will enter a low battery protection state. At this time, if you restart the operation, the fascia gun will not work properly.)

Solution: Keep the device turned off, connect the charger, make sure the battery is fully charged, and then try using it again.

*** Operational troubleshooting

1. Product not working: Check if the battery is sufficient, it is recommended to use it after fully charging it.
2. Slight abnormal noise during use: The new machine running in period s a normal phenomenon. If there is a very obvious abnormal noise, please send it for inspection.
3. The battery cannot be charged: Battery aging or poor contact, please contact customer service or sales personnel

*** Statement:

To reduce the risk of malfunction, electric shock, injury, fire, death, and equipment damage, please read all instructions and warning prompts on the product in this manual before using this product.

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