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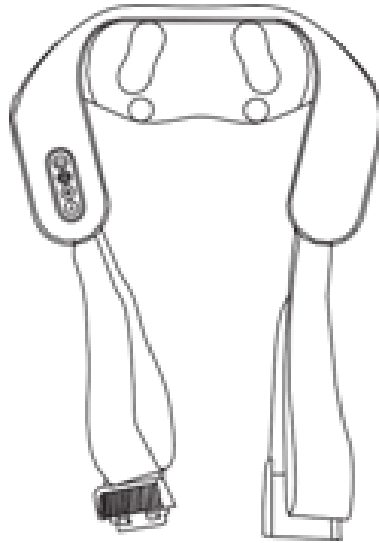
# INSTRUCTION MANUAL


NECK AND SHOULDER MASSAGER

MODEL: MK366

Read this instruction manual carefully before using

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## **Product Instructions**
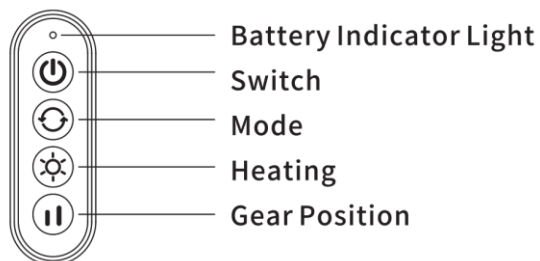
## **Product Parameters**

| Product name               | Rated voltage | Rated power | Product size        | Product weight |
|----------------------------|---------------|-------------|---------------------|----------------|
| Neck and Shoulder Massager | DC 5V         | 25 W        | ≈ 23 x 13.5 x 16 Cm | ≈ 1.18 Kg      |

## **Operating Instructions**

The product is powered by two 18650-7.4V lithium batteries, which are charged using USB 5V-2A, Type-C. When the product's battery level is low, the red indicator light flashes 6 times and then automatically shuts down.

When charging, the product cannot be turned on for use, and the indicator light is in a red breathing state. When fully charged, the indicator light turns white and flashes. The product is scheduled for 10 minutes.

### ***\*\*\*Correct usage method***

#### **\* Preparation before use**

Please check if the product and accessories are complete and undamaged.

Press and hold the on button for 2 seconds, and the product starts working.

#### **\*Mode selection**

The default mode for startup is mode 1. Short press the mode key to switch between mode 1 and 2.

#### **\*Hot compress mode**

The default hot compress mode is high when starting up. Short press the hot compress button, and then go to low compress mode- close hot compress mode to cycle.

#### **\*Strength mode**

The default gear is first gear when starting up. Short press the gear button to switch between the first and the second gear.

#### **\*Product shutdown**

Automatic shutdown: After a total of 10 minutes of use, the product stops working.

Manual shutdown: Press and hold the on button for 2 seconds to stop the product from working.

## **Precautions For Product Disposal**

The lithium-ion battery used in this product is a valuable resource that can be recycled and regenerated. Please take out the lithium-ion battery when it is used and discarded and hand it over to a professional organization that can recycle rechargeable Battery recycling for treatment.

## **Trouble shooting**

| Fault Phenomenon       | Possible Causes                                                                | Handling Method                             |
|------------------------|--------------------------------------------------------------------------------|---------------------------------------------|
| Unable To Start        | The battery is dead                                                            | Turn on the product after fully charging it |
| Stopping<br>During Use | The time for reaching the scheduled time has been used                         | Please turn off the power and restart       |
|                        | Excessive product load (such as forcefully stepping or stacking heavy objects) | Please reduce the load and use it           |
| Sound During Use       | Emitted from the mechanical structure of the product itself                    | Normal state                                |

## **Emergency Handling Methods**

1. If there is any damage to the power cord of this product, please stop using it immediately to avoid danger.
2. If the malfunction is caused by improper use, please do not disassemble it yourself. Please unplug the power plug and contact customer service for repair.
3. If the wires or plugs of the machine are damaged, please stop using them immediately.
4. If smoke or fire occurs during use, please unplug the plug from the socket as soon as possible and notify customer service.
5. If the leather or fabric surface is damaged and the machine parts are exposed, please stop using immediately and contact customer service for repair to avoid danger.

## **Maintenance**

For daily maintenance, please use neutral detergent and soft cloth to scrub.

Do not use engine oil, organic solvent or other chemicals to clean. Do not rinse in water.

Please gently press the switch during operation to avoid excessive force.

Avoid heavy pressure and sharp objects from scratching the surface of the product.

After use, please keep it properly and place it in the original box or in a dry and low-temperature place.

If any irreparable faults are found during use, please immediately cut off the power and stop using, and contact customer service. Please do not disassemble by yourself.

Please unplug the power promptly when charging is complete.

## **Note**

Do not pull the power cord to avoid causing wire breakage or poor contact; When unplugging the power plug, hold the plug and unplug it. In case of abnormal situations such as the entire machine not working, immediately stop using and turn off the power.

Do not allow infants or children to use or play on the product.

When feeling unwell during use, please discontinue use.

Individuals who are pregnant or have menstrual periods, high fever, malignant tumors, heart disease, incomplete surgery, acute illness, or other illnesses should not use it according to medical advice.

●This product is prohibited from being used in the following places: bathrooms, kitchens, and other places with high temperature, high pressure, humidity, and near water and fire sources.

Please stay away from sharp objects. When the outer cover is damaged, please do not use it.

Please contact customer service for repair.

It is strictly prohibited to use the machine leather cover when it is broken to prevent hair from getting entangled in the machine and causing serious personal injury.

Please shut down or unplug the power cord from the power outlet after use or before cleaning.


If the power cord is damaged, please contact customer service for repair and handling in a timely manner to avoid danger.

It is necessary to use the built-in charger of this product for charging, otherwise the product maybe damaged or even ignited.

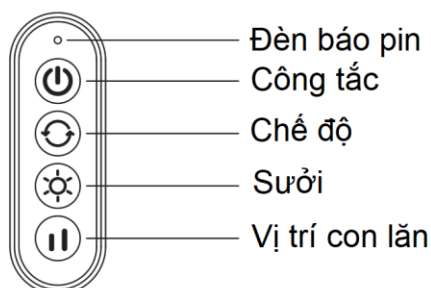The device has a heating surface and must be used with caution by people who are not sensitive to heat.

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