LocknLock

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3.7V, 8W, 250ML

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LocknLock

INSTRUCTION MANUAL

RECHARGEABLE ORAL IRRIGATOR ENR154WHT
3.7V, 8W, 250ML

Read instruction manual before using to use product in safe and appropriate way.
Keep this manual in the proper place in order that we can take it when we need it.

CONTENTS

Product specification	12
Safety instruction	13
Product structure	14
Before first use	15
Operation	15
Care and maintenance	17
Troubleshoot	19

PRODUCT SPECIFICATION

Product name	Rechargeable oral irrigator
Item code	ENR154WHT
Voltage	3.7V
Power	8W
Input	5V  1A
Capacity	250ml

SAFETY INSTRUCTION

All information in this section is provided to instruct you in the correct and safe operation of this device so that you can avoid injury to yourself and others as well as damage to property. Read all instructions before using.

DANGERS


1. Toddlers, infants, and people who are unable to operate should not be allowed to use this appliance.
2. Please keep it out of reach of children.
3. Do not try to disassemble or repair it.
4. Do not place or store the product where it can fall or be pulled into tub or sink.
5. Do not put the product in fire or heat it; charge, use or store it in high temperature environment, as it contains a built-in rechargeable battery.
6. Do not use the product while bathing.
7. Do not charge when the body is wet.

WARNINGS

1. This appliance can be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
2. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
3. The battery must be removed from the appliance before it is scrapped.
4. Do not fill the tank with water warmer than 40°C, otherwise it might lead to burns.
5. If the charging cable is damaged, it shall be replaced by customer service center or a similar qualified person in order to avoid hazard.
6. Please do not use the appliance if nozzle is lost.
7. Do not submerge the device in water.
8. Do not use it for any purpose other than oral cleaning.
9. Keep packing bags and nozzle of this product away from children. The nozzle is small enough to be swallowed and can be lodged in the throat.
10. People who experience tooth or gum pain may not be able to use this device and should consult the dentist.
11. If the gums bleed profusely after use of this product or still after 2 weeks of use bleeding occurs. Please stop using this product and consult your dentist.
12. Do not press the nozzle too firmly against your teeth or gums. Doing so may injure your gums.
13. Do not share nozzles with family members or other personnel. It may cause infection or inflammation.
14. Certain types of mouthwash may damage the equipment by causing the housing or water tank cracked. Therefore, do not use mouthwash in the water tank.
15. After use, clean the product and nozzle and allow it to air dry. Otherwise, it may cause damage to health due to the growth of mold and bacteria.
16. Stop using the product immediately and contact customer service center in these below situations:
 - The product becomes abnormally hot.
 - The product emits a burning smell.
 - The product makes abnormal sounds during use or charging.

PRODUCT STRUCTURE

PRODUCT DESCRIPTION

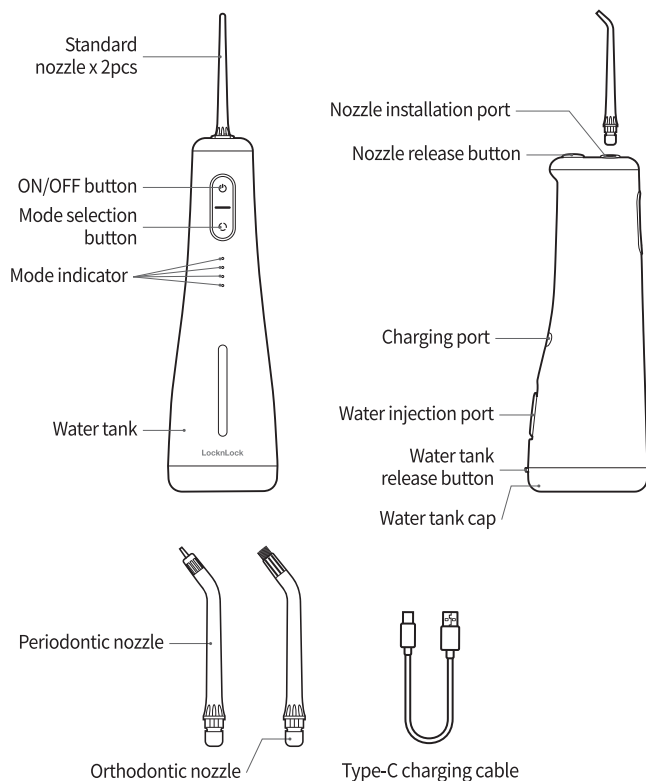NOTICE

There may be some liquid remaining inside the main appliance before use (tank, water suction tube, etc.), this is remaining distilled water that was used for the performance (waterflow) test, so there are no health concerns.

NOZZLES

1. Standard nozzle x 2pcs: effectively remove food residue, tartar, plaque in the teeth, and prevent tooth decay and dental calculus.
2. Periodontal nozzle x 1pc: effectively remove bacteria from periodontal pockets, massage the gums, prevent gingivitis and gingival atrophy.
3. Orthodontic nozzle x 1pc: ideal for cleaning around dental bridges, braces and crowns.

BEFORE FIRST USE

1. Please read this manual before operating the appliance and keep it in a safe place for future reference.
2. Remove all package materials before first use.

OPERATION

HOW TO USE

1. Attach nozzle onto the unit. (Fig. 1)
2. Open the water injection cap, put the water injection port horizontally and fill the water tank, then fasten the water injection cap. (Fig. 2)

TIPS

1. It's recommended to use distilled or purified water.
2. Water filled in the water tank should not exceed 40°C.
3. Press the MODE button to select the mode you desired. During irrigating, you can change the operation mode.

Clean: daily clean use

Soft: gentle cleaning for sensitive gums

Pulse: deep cleaning and massaging the gums

DIY: modify the desired water pressure, kindly follow the below instruction:

- a. Press ON/OFF button.
- b. Choose DIY mode, then keep pressing the MODE button.
The water pressure is increased 7 levels from the softest to the strongest.
- c. When release the MODE button, the related water pressure will be selected.

4. Hold the unit in an upright position with nozzle tip point toward your teeth.
5. Keep your mouth opened slightly so that the water can flow out.
6. Press ON/OFF button to start irrigating. (Fig. 3)
 - a. Directing the water jet flow at a 90-degree angle to the teeth and gums.
 - b. Directing the water jet slowly along the teeth.
 - c. Directing the water jet along the gum line with the nozzle tip at a 90-degree angle to the gum tissue.
 - d. Cleaning between the teeth and around any braces, crowns or bridgework thoroughly.

7. Press ON/OFF button and then press the nozzle release button to replace other nozzles when you need. (Fig.4)

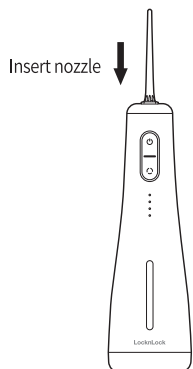
Fig. 1

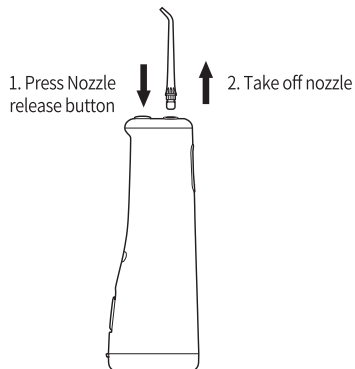
Fig. 4

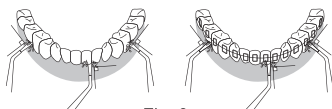
Fig. 3

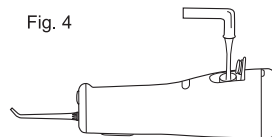
Fig. 2

NOTICE

When used for the first time or not used for a long time, after the product is turned on, it might take a few moments for the nozzle to spray water. At this time, please remove the nozzle first, and then add an appropriate amount of water to the nozzle installation port, and then install the nozzle before using this product again.

HOW TO CHARGE

The number of indicators is on after the product is shut down showing the battery status. When the four indicators blink at the same time after the product is turned off, it indicates that the battery is low, please follow the steps below to charge:

1. Open the charging port cover and insert the appliance plug securely into the charging port (make sure there is no liquid remaining inside the charging port before inserting the appliance plug).
2. Insert the charging adaptor securely into the power outlet.

3. During charging, indicator light will flash. When all the 4 indicators are on, it indicates that the product is fully-charged. It takes about 4 hours for this product to be fully-charged.
4. Install the charging port cover again after charging finished.

NOTICE

1. For the first time charging or not used for a long time, the indicators might not be on in the first 5 minutes, which is caused by insufficient battery. This is a normal phenomenon.
2. The product needs to be charged in the following cases
 - First use
 - The product is not used for more than 3 months.
 - When using, it's found that the water pressure is low.
 - When there are low power indicators.
3. The product might become hot during charging, which is a normal phenomenon.
4. The recommended charging ambient temperature is from 5°C to 35°C. When the room temperature is over or below than the recommended, the battery performance may be degraded.

CARE AND MAINTENANCE

1. Use only water or neutral cleaning agent to clean the device. Never use caustic or abrasive agents (e.g. vinegar or lime scale removers), since they may damage the device.
2. Do not use hot water of over 50°C for cleaning.
3. If the product is not used for a long time, the product should be stored in a cool, dry place away from direct sunlight.
4. Do not use the dryer or the heater to dry the product, otherwise it may cause component failure or deformation.

AFTER USING THE APPLIANCE (Fig.5)

1. Open the bottom water tank cap and pour out the remaining water.
2. Turn the product on for a period of time to drain water inside.
3. Press the ON/OFF button to shut down the product, and then install the water tank bottom cap again
4. Wipe the product with a dry, soft cloth.

NOTICE

To prevent bacteria propagating and prolong product lifetime, make sure there is no remaining water inside the unit. If this appliance is not used for long time, be sure to wipe and dry out the appliance completely before storing it.

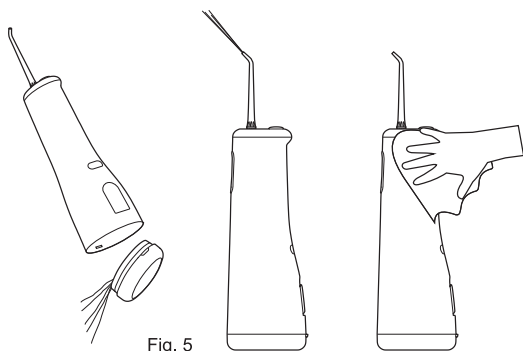
Fig. 5

MAIN UNIT

1. Remove the dirt by wiping using a cloth with liquid soap and wipe off using a dry cloth.
2. Do not submerge the main unit in water to clean.
3. Dry after washing with water.

WATER TANK CAP

1. Wash with water and wipe with soft cloth.
2. Dry the water tank.

NOZZLES

1. Wash with water and wipe with soft cloth.
2. It is recommended to replace the nozzle every 6 months.

DISPOSAL

The rechargeable battery inside your appliance is easy to remove for recycling. At the end of the product's life and prior to disposal, please remove the rechargeable battery from the handle by following the instructions below.

During below operation, be sure to protect your eyes, hands, fingers, and the surface on which you work.

1. Remove the bottom water tank cap of the water tank (① in Fig. 6) and the charging port cover (② in Fig. 6).
Open the top cover with a pointed hard object (③ in Fig. 6) and remove the nozzle release button (④ in Fig. 6), spring (⑤ in Fig. 6) and nozzle lock (⑥ in Fig. 6)
2. Use a screwdriver to remove the two screws at the top of the product (Fig. 7)
3. Remove the product cover from the top (Fig. 8)
4. Cut the battery connection wire by scissors and take out the battery (Fig. 9)

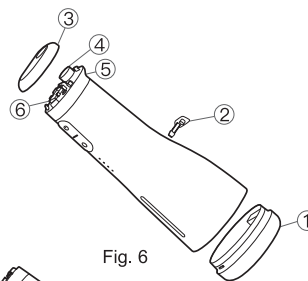
Fig. 6

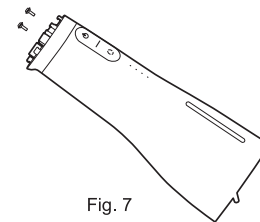
Fig. 7

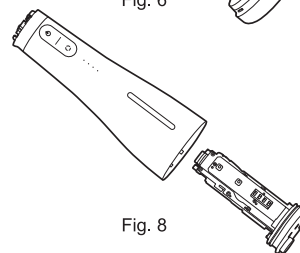
Fig. 8

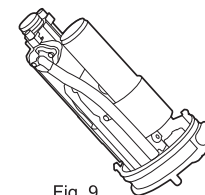
Fig. 9

TROUBLESHOOTING

Problem	Possible cause	Solution
It does not work	The appliance has just been purchased or it has been left unused for more than 3 months	Charge the appliance for at least 4 hours
It only works for a few minutes even after charging	The battery has reached the end of its life cycle	Contact customer service center for checking
	Charging time was insufficient	Charge the appliance for at least 4 hours
Water pressure is very low	There is dirt stuck in the nozzle	Replace the nozzle
	The nozzle is deformed	Replace the nozzle
Water does not come out	The water tank is empty	Fill up the water tank
	The appliance is tilted too much while in use	Use the appliance in vertical state