CHOICE Dr.

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CHOICE Dr. **DC-004**

CHOICE Dr.

Thank you for choosing and supporting CHOICE Dr.!

We hope that you are satisfied with our product and service, and that we have the opportunity to serve and accompany you and your family.

If you have any feedback or reviews regarding the quality of our product and service, please contact our hotline at **1800 8379** so that we can receive and promptly address them.

In addition, please scan the **Zalo OA** code to activate the electronic warranty, receive additional gifts, and enjoy promotion and **24/7** support from CHOICE Dr.


Thank you once again and we are happy to have take care of you and your family's health!

Choice Dr.

PRODUCT SPECIFICATIONS

Product Name: Leg Massager	Model: DC-003
Power: 8W	Operating Voltage: 3.7V (maximum 4.2V)
Charging Voltage: 5V - 2A	Product Weight: 340 gram
Product Dimensions: 460x540x260 mm	

PRODUCT STRUCTURE

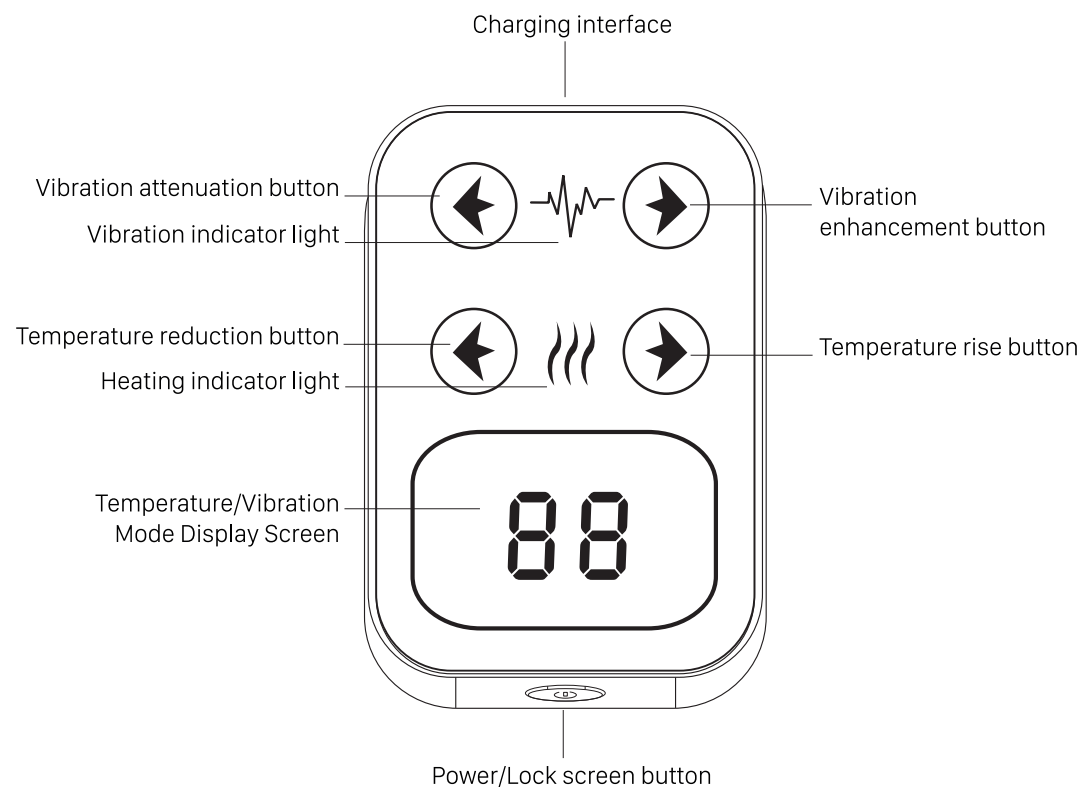OPERATING INSTRUCTIONS

01

Press and hold the power button to start the device.
A beep will sound, the indicator light will turn white, and the default mode will start at the first low-temperature level.
(If the device cannot be turned on, please charge it and try again.)

02

- Press the right heating button to increase the temperature or the left button to decrease it.
- When the maximum temperature is reached, press the right button again to maintain the temperature, and a beep will sound.
- If the heating function is off, press the left button to maintain the no-heating mode — the buzzer will beep twice.
- Temperature levels displayed on screen:
45°C = Level 1 50°C = Level 2
55°C = Level 3 60°C = Level 4
65°C = Level 5 00°C = No heating

(Testing environment and conditions comply with GB4706.1-2005 and GB4706.8-2008 standards.)

03

- Press the right vibration button to increase vibration intensity or the left button to decrease it.
- When reaching level 3, press the right button again to keep the vibration constant — a beep will sound.
- When the vibration is turned off, press the left button to keep the no-vibration mode — the buzzer will beep twice.
- Vibration modes displayed on screen:
01 = Mode 1 02 = Mode 2
03 = Mode 3 00 = No vibration

04

During operation, to prevent accidental touches, you can short press the power button to lock or unlock the screen.
When locked, the display will turn off (only the battery indicator remains white).

05

Press and hold the power button again.
The buzzer will beep twice, and the device will power off.

CLEANING & MAINTENANCE

- Ensure the device is turned off and unplugged before cleaning.
- Wipe the device gently with a dry cloth.
- Keep the product clean and store it in a dry place.
- Remove dust using a soft cloth dampened with mild soapy water, then wring the cloth before wiping again.
- Do not use: alcohol, paint thinner, gasoline, or nail polish remover, as these may damage the surface or cause discoloration and cracking of components.

SAFETY WARNINGS

- Massage should feel comfortable and relaxing. Stop using immediately and consult a doctor if you feel pain.
- Do not use non-standard power cables, as they may cause overheating or fire.
- Do not use the device on injured or damaged areas.
- Do not use after taking painkillers, sedatives, or alcohol.
- Avoid using on numb or insensitive skin or in poorly ventilated environments.
- Do not use within 30 minutes before or after meals.